




EDUCATION, TRAINING, AND 24/7 SUPPORT

Our dedicated support team provides a wealth of knowledge, resources and technical training opportunities, including hands-on product training for clinical personnel, an interactive training app, online videos, and 24/7 support.

More than just tubing, the ambIT[®] cassette is the pumping mechanism behind the ambIT[®] infusion system that attaches to the pump driver. Each single-use cassette, with its built-in rotary peristaltic pump technology, is manufactured to exact tolerances, providing accuracy and ease of use.

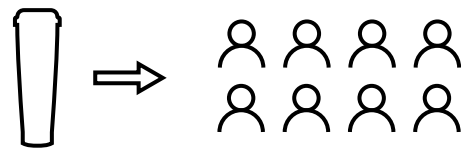

There are inherent risks in all medical devices. Please refer to the product labeling for Indications, Cautions, Warnings, and Contraindications. Failure to follow the product labeling could directly impact patient safety. Physician is responsible for prescribing and administering medications per instructions provided by the drug manufacturer.

For more information please visit
avanospainmanagement.com
1-800-448-3569
844-4-AVANOS
painmanagement@avanos.com

For more information on the full line of
Avanos pain management products, visit
avanospainmanagement.com.
© 2020 AVNS All rights reserved.
COPY-04632 08/20

AVANOS



AMBIT[®] FAMILY OF PUMPS
An ambulatory infusion pump for every need

AVANOS

avanospainmanagement.com



ambIT[®] pumps are simple to set-up, light-weight, programmable, and have a volumetric accuracy of +/-6%.

Suitable for multiple patient care areas, procedures and treatments, ambIT[®] pumps provide a variety of infusion rate options.

ambIT[®] pumps are simply powered by two "AA" batteries and can be programmed in minutes.

Features include:

- Portability
- Durability
- Reusability
- Light-weight
- Accuracy
- Easy to Program



Accurately addresses pain.



Programmable settings enable physicians to achieve a patient-specific dose allowing for a targeted pain management experience.



ambIT[®] bolus button allows patients to self-administer small doses of medication to help reduce break-through pain and improve patients' comfort.



Volumetric accuracy $\pm 6\%$ ensures patients receive the appropriate amount of medication.



Physicians can reassess and adjust pain management plans as needed.



To learn more about ambIT[®] call
1-800-448-3569
844-4-AVANOS
 or email painmanagement@avanos.com



PRESET™

A disposable device with a fast and easy way to program. Provides four pain management protocols to choose from and a fifth customizable option.

ONLY AVAILABLE IN DISPOSABLE FORMAT

PN: 220426



PIB•PCA

Combines the PCA pump functions with continuous and intermittent therapies, allowing flexibility based on patient need.

PN: 220244



PCA

Provides a simple and accurate solution for post-operative local anesthetic infusions.

PCA has three profiles: basal only, basal plus patient-controlled bolus and basal plus bolus.

PN: 220262



PIB•PCA

Combines the PCA pump functions with continuous and intermittent therapies, allowing flexibility based on patient need.

PN: 220501



KIDS•PIB

Provides the same functions of the PIB•PCA, plus an additional feature. With the KIDS•PIB, a bolus is administered through an external remote that detaches from the pump when bolus is not being administered.

PN: 220540



PIB•IL

Uses the same features and flexibility you find with the PIB•PCA, with a programmable interval limit feature.

PN: 220527