

MEET PATIENT'S NEEDS BY COMBINING REGULARLY OCCURRING DOSES AND MANUAL BOLUS OPTIONS

PIB+PCA Mode

Regularly occurring doses of medication while allowing patients to request a manual bolus in between intermittent doses.

The ambIT® PIB•PCA Pump is a safe, simple, multipurpose infusion therapy system

Effective Safety Standards:

- Programming lockout
- Bolus lockout
- Anti-free flow
- Program review
- Device internal self-check
- Latex-free and DEHP-free fluid path

Education, Training, + 24/7 Support

Our dedicated support team provides a wealth of knowledge, resources and technical training opportunities, including hands-on product training for clinical personnel, an interactive training app, videos, and 24/7 support.



Three Unique Programs:

- PCA** Patient Controlled Analgesia
Continuous basal rate and optional patient-controlled bolus
- PIB** Programmable Intermittent Bolus
Regularly occurring doses of medication
- PIB + PCA** Combined PCA and PIB
Regularly occurring doses of medication and manual bolus options

For service or to place an order:

1-800-448-3569


844-4-AVANOS

painmanagement@avanos.com

For more information on the full line of Avanos pain management products, visit avanospainmanagement.com.

Volumetric Accuracy	+/- 6%	
Dose Volumes	0 to 50 ml in 0.1 ml steps	
Dose Interval Times (hr:min)	00:01 to 24:00 in 1 minute increments	
Basal Flow Rates	0 to 20 ml/hr in 0.1 ml/hr steps	
Stroke Volume	50 microliters	
Bolus Delivery Rate*	210 ml/hr	
Bolus Dose Volume	0 to 50 ml in 0.1 ml steps	
Bolus Lockout Times (hr:min)	00:01 to 24:00 in 1 minute increments	
Volumes to be Infused (ml)	0 to 1000	
Pump Mechanism	Microprocessor Controlled Rotary-Peristaltic	
Dimensions	2.16 in. x 1.4 in. x 6.875 in. (55 mm x 36 mm x 175 mm)	
Weight	6.4 ounces (181.4 grams) with batteries	
Power Supply	2 AA 1.5V batteries	
Alarms & Alerts	<ul style="list-style-type: none"> • Run Indicator Light • Bolus Infusing • Pause Indicator • Occlusion Downstream • Detached Cassette 	<ul style="list-style-type: none"> • Low Battery • Dead Battery • Malfunction • Boluses Requested • Infusion Complete
Reports	<ul style="list-style-type: none"> • Volume Infused • Boluses Delivered 	<ul style="list-style-type: none"> • Boluses Requested • Elapsed Time
Delivery Profiles	<ul style="list-style-type: none"> • PIB+PCA • PCA • PIB 	

*rate varies

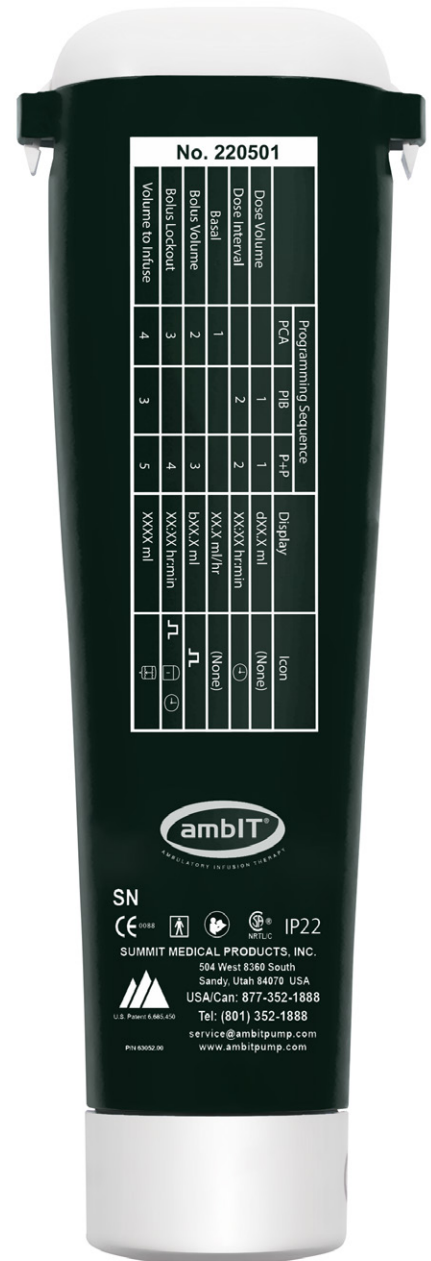

The ambIT® PIB•PCA Pump follows internationally recognized volumetric accuracy standards (IEC 60601-2-24).

There are inherent risks in all medical devices. Please refer to the product labeling for **Indications, Cautions, Warnings, and Contraindications**. Failure to follow the product labeling could directly impact patient safety. Physician is responsible for prescribing and administering medications per instructions provided by the drug manufacturer.