

# A SAFE, SIMPLE, MULTIPURPOSE INFUSION THERAPY SYSTEM

ambIT® zeros in on the most essential features to ensure patient's needs are met in a **simple, safe, and adaptable** manner.

#### ambIT® PCA:

A patient-controlled, electronic bolus pump that maximizes functionality and versatility to meet patient infusion needs.

#### ambIT® PreSet:

Fast and easy programming with four pain management protocols to choose from and a fifth option to customize your own.



#### Flexibility and Convenience:

- Small and lightweight design
- Easily concealed in custom-built carry pouch with view-window
- Easy-press bolus button
- Titratable infusion rates
- History reports
- Disposable cassette

#### Effective Safety Standards:

- Programming lockout
- Bolus lockout
- Anti-free flow
- Program review
- Device internal self-check
- Latex-free and DEHP-free fluid path


For service or to place an order:

**1-800-448-3569**

**844-4-AVANOS**

**[painmanagement@avanos.com](mailto:painmanagement@avanos.com)**

For more information on the full line of Avanos pain management products, visit [avanospainmanagement.com](http://avanospainmanagement.com).

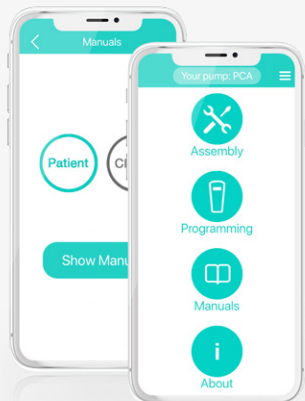

### Education, Training, + 24/7 Support

Our dedicated support team provides a wealth of knowledge, resources and technical training opportunities, including hands-on product training for clinical personnel, an interactive training app, videos, and 24/7 support.

Volumetric Accuracy	+/- 6%
Basal Flow Rates	0.1 to 20 ml/hr (increments vary by model)
Stroke Volume	50 microliters
Bolus Delivery Rate	100 ml/hr
Bolus Dose Volume (ml)	0.1 to 20 ml
Bolus Lockout Times (hr:min)	00:01 to 12:00 (increments vary by model)
Volumes to be Infused (ml)	25 to 1000 ml
Pump Mechanism	Microprocessor Controlled Rotary-Peristaltic
Dimensions	2.16 in. x 1.4 in. x 6.875 in. (55 mm x 36 mm x 175 mm)
Weight	6.4 ounces (181.4 grams) with batteries
Power Supply	2 AA 1.5V batteries
Alarms & Alerts	<ul style="list-style-type: none"> <li>• Run Indicator Light</li> <li>• Bolus Infusing</li> <li>• Pause Indicator</li> <li>• Occlusion Downstream</li> <li>• Detached Cassette</li> <li>• Low Battery</li> <li>• Dead Battery</li> <li>• Malfunction</li> <li>• Boluses Requested</li> <li>• Infusion Complete</li> </ul>
Reports	<ul style="list-style-type: none"> <li>• Volume Infused</li> <li>• Boluses Delivered</li> <li>• Boluses Requested</li> <li>• Elapsed Time</li> </ul>
Delivery Profiles	<ul style="list-style-type: none"> <li>• Basal Rate + Bolus</li> </ul>

The ambIT<sup>®</sup> PIB•PCA Pump follows internationally recognized volumetric accuracy standards (IEC 60601-2-24).

There are inherent risks in all medical devices. Please refer to the product labeling for **Indications, Cautions, Warnings, and Contraindications**. Failure to follow the product labeling could directly impact patient safety. Physician is responsible for prescribing and administering medications per instructions provided by the drug manufacturer.